

Crown - Master[®]

Order Number I12019

Diaduran[®], the bend resistant, heat resistant fibre-glass plier.



For a tighter fit, simply reduce the size of the knob.




For the exact smoothing, finishing and polishing of the margins.



Use this hold for the precise finishing of labial or incisal surfaces.




Work on the entire crown by simply turning the plier.

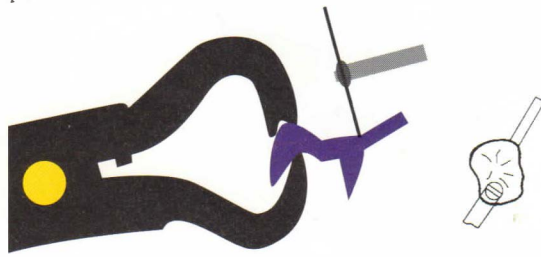


The ideal tool for holding bridges or crowns while separating sprues with a disc.



Makes it possible to work on the entire occlusal surface.

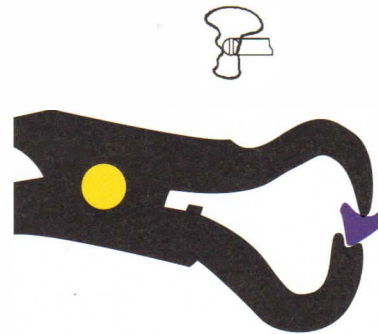

Easy positioning of inlays.

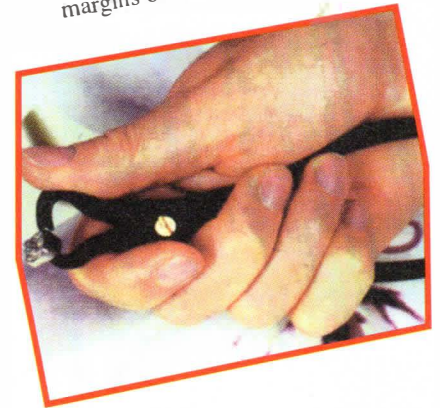

Provides relative ease for smoothing and polishing the margins of molar crowns.



The groove in the handle of the plier makes it possible to work on a single crown. Place crown in the groove provided for finishing, **dust off**, polish.

Pfingst
AND CO., INC. SINCE 1905